

Ingenuity in
infection control
since 1830



Little Sister SES 2010 Non-vacuum Autoclave

User Manual

Thank you for choosing Eschmann

Eschmann design, manufacture and supply a range of market leading products including benchtop autoclaves, accessories, washer disinfectors and surgical suction.

For further information visit our website: www.eschmann.co.uk.

Please read the manual before installing or using the product for the first time. Keep the manual handy for quick reference. Always make sure that the manual is available for the next user/owner of the product.

The product(s) described by this manual must only be operated and used by qualified personnel trained in the use of this equipment. Contact Eschmann with your training requirements.

Pay particular attention to the safety notes, cautions and warnings provided in the text, and also to those displayed on the product labels.

This product must be used, installed and maintained in accordance with the procedures given in this manual. Failure to do so could result in injury to patients and/or users, or damage to the equipment.

Eschmann products must be serviced by Eschmann trained engineers only. Failure to do so may invalidate the warranty.

Contents

1	Safety Warnings and Usability	3	12	Menus	17
1.1	Safety checks/ tasks	3	12.1	Cycle Selection Menu	17
1.2	Potential hazards	3	12.2	Settings Menu (Managers only)	19
1.3	Limitations on use	3	12.3	Records Menu	22
1.4	Electrical safety	3	13	User security	24
1.5	General safety	3	13.1	Overview	24
1.6	Warnings, cautions & notes	3	13.2	Autoclave manager	24
1.7	Usability	4	13.3	Setting a PIN code	25
Fig 1	Main Components	5	13.4	Selecting the security level	26
Fig 2	Typical Start Screen	5	14	Cycle Logger application	27
2	Introduction and Installation	6	14.1	Install the software	27
2.1	Introduction	6	14.2	Create a PC back-up	27
2.2	Installation	6	14.3	Viewing files	27
2.3	Unpacking	6	14.4	Printing files	28
2.4	Autoclave positioning	6	14.5	Connecting a printer	28
2.5	Mains connection	6	Appendix 1 : Daily/weekly testing	30	
2.6	Information/ safety symbols	6	Daily Testing	30	
3	Controls	7	HTM01-01 (HTM2010/DB2002(06)) testing	30	
3.1	Touch screen & menu operation	7	Appendix 2 : Loading the autoclave	31	
3.2	The Keypad	7	Unwrapped solid items	31	
3.3	The Progress Screen	8	Using the tray lifter	31	
3.4	Notification screens	8	Load capacity	31	
3.5	Autoclave manager	8	Typical loads	32	
3.6	Water Quality	8	Appendix 3 : Maintenance	33	
4	Preparing the autoclave	9	Cleaning and care	33	
5	Loading the autoclave	9	Autoclave fuses	33	
6	Starting and running a cycle	10	Mains cable fuse	33	
7	Cycle completion and unloading	11	Waste electrical & electronic equipment (WEEE) regulations	33	
8	How to stop the cycle	12	Service, Calibration & Certification	33	
9	After use	13	Appendix 4: Troubleshooting	34	
9.1	Emptying the water reservoir	13	Trouble shooting table	34	
9.2	Switch off	13	Advisory code table	34	
9.3	Cleaning	13	Appendix 5: Technical data	37	
10	Service Screens	14	Environmental Requirements	37	
10.1	Replace the door seal	14	BS EN 13060 Type Tests	39	
10.2	Service required	14	Load Chart	39	
11	Advisory Screens	15	Appendix 6: Accessories	41	

1 Safety Warnings and Usability

WARNING!

Only use the autoclave, as specified in these instructions. Eschmann are not responsible for a malfunction or reduced level of protection provided by the equipment when not used as specified. Never tamper with, bypass or interfere with any of the safety features.

Only use Eschmann accessories and mains leads. Any spare parts required must be supplied by Eschmann. They must be fitted in accordance with the Service Manual or instructions supplied by Eschmann, or Eschmann trained engineers

Be careful not to splash water on the top of the autoclave. If water is poured onto the autoclave, or a leak is found, disconnect from the mains socket, dry thoroughly and allow time to dry

Servicing and repairs **MUST** be performed by Eschmann, or Eschmann trained engineers

1.1 General safety

Daily

- Drain the water reservoir at the end of each day.
- Clean the door seal and chamber face every day with Eschmann recommended wipes.
- Check for escape of steam or water during a cycle (other than from the water reservoir lid)
- Replace the water reservoir lid when filling with water is complete. Refer to Section 3.6, Water Quality.

Weekly

- Check the door seal and chamber face for any signs of damage.
- Check condition of the mains lead, mains plug and outer covers for signs of damage.

1.2 Potential hazards

- Autoclaves operate with steam at high pressure and temperature. Locations of potential hazards are marked on the autoclave. Refer to Section 2.6.
- Take care to avoid contact with any residual steam or hot water as the autoclave door is opened.
- Take care to avoid contact with the hot internal surfaces of the chamber and door. Use the Eschmann tray lifter to load and remove trays from the chamber.
- Water in the reservoir can be very hot. Take care when draining to avoid contact with the hot water.
- The water reservoir lid, on top of the autoclave, vents hot steam during normal operation. Take care to avoid contact. Do not cover the water reservoir lid or remove it during operation.
- If a warning screen is displayed, the load **MUST** be treated as non-sterile. Refer to Appendix 4, Troubleshooting. In this case, **DO NOT** switch off power until discharge of steam into the reservoir has stopped. Do not attempt to open the chamber door until the unit has cooled.

1.3 Limitations of use

- **THIS IS AN N TYPE AUTOCLAVE, SUITABLE ONLY FOR UNWRAPPED, SOLID ITEMS.**
- **DO NOT** process hollow, porous, pouched or wrapped items.
- **DO NOT** process liquids.
- **Never** use trays or cassettes without perforations.
- **DO NOT** use near flammable materials, or gases.
- **Items MUST be cleaned before sterilisation in the autoclave. If instruments are not cleaned first, this may compromise the sterilisation process.**
- **Sterilised items are intended for immediate use.**

1.4 Electrical safety

- This equipment must be earthed.
- Eschmann recommend the use of a suitable RCD.
- Always switch the autoclave off and disconnect from the mains electrical supply before renewing fuses, checking and cleaning the autoclave.
- No attempt must be made to service the autoclave internally.

1.5 General safety

- Ensure the water reservoir lid is fitted to the top of the autoclave, prior to running.
- **DO NOT** fill the water reservoir during operation.
- **DO NOT** cover the water reservoir lid.
- **DO NOT** use the autoclave if the door seal, or outer covers are damaged, loose or missing.
- **DO NOT** place heavy items or containers with liquids (e.g. cups etc.) on top of the autoclave.
- **DO NOT** cover the ventilation grills located on the back and underside of the autoclave.
- **DO NOT** press the touch screen too hard, or use sharp objects to press the screen.
- **DO NOT** use flammable liquids, abrasive powders, chemicals, or solutions containing chlorine to clean the autoclave.
- Potentially contaminated waste materials produced during cleaning and disinfection should be handled in accordance with local procedures and national legislation for the disposal of potentially contaminated waste.

1.6 Warnings, cautions & notes

Various warnings, cautions and notes are made throughout this manual. Each of these carries a special meaning and should be read carefully.



A **WARNING** is given when the safety of the patient or user may be involved. Disregarding this information could result in injury to the patient or user.

A **CAUTION** is given when special instructions must be followed. Disregarding this information could result in permanent damage to the product.

A **NOTE** provides specific information that makes important instructions clear.

1.7 Usability


Intended User	Trained medical professionals. The product is not used by the patient.
Medical Purpose	Sterilisation of medical instruments using saturated steam.
Contraindications	<p>There are no direct contraindications. The autoclave does not provide direct treatment to a patient. Indirect contraindications include cycle failure, incorrect loading and incorrect servicing/ repair, leading to a non sterile load.</p> <p>The risk of non-sterilisation has been assessed through risk management and mitigated through product design, verification/ validation product testing, and instructions for use, in accordance with the standards listed in Appendix 5.</p> <p>In the event of cycle failure, the autoclave is designed to alert the operator that the cycle has failed and the load is not sterilised.</p> <p>To ensure a sterilised load, the autoclave must be operated by trained medical staff and used in accordance with the instructions and safety warnings contained within this manual.</p> <p>Servicing and repairs must be performed by Eschmann, or Eschmann trained engineers.</p>
Patient Population	The autoclave sterilises medical instruments and therefore does not provide treatment directly to a patient.
Intended use	<p>Designed and tested for compliance to BS EN 13060, the product is a bench top small steam steriliser for the sterilisation of medical instruments using saturated steam.</p> <p>The autoclave has an 11 litre chamber and operates using an N Type cycle at a nominal temperature of 134°C, either with or without drying.</p> <p>The product is intended for the sterilisation of unwrapped solid instruments.</p> <p>Refer to Appendix 2 for full details of the loading requirements and capacity.</p> <p>Refer to Appendix 5 for Technical data.</p> <p>The product is intended to sterilise clean items only.</p> <p>The product is not suitable for hollow, porous, pouched or wrapped items.</p> <p>The product is not suitable for sterilising liquids or pharmaceutical products.</p> <p>The product is not suitable for sterilising items which are not rated for steam sterilisation at the stated cycle temperatures.</p> <p>The product has a technical lifetime of 7 years.</p>
Device Application	The autoclave is suitable for the sterilisation of solid, unwrapped medical instruments in the following applications: Hospitals, doctors surgeries, medical centres, dental practices, podiatry practices, veterinary practices, hair salons, and tattooists.
Device Classification	Devices specifically to be used for sterilising medical devices are classified as Class IIb as derived from Rule 15 of Annex IX of the EU Medical Devices Directive 93/42/EEC and the UK Medical Devices Regulations 2002 (SI 618), as amended by the EU Exit Regulations 2019 (SI 791) and 2020 (SI 1478).
Frequently used functions	<ul style="list-style-type: none"> • Switching the autoclave on/ off • Filling, emptying and changing the water in the reservoir. • Opening and closing the door • Loading and unloading the autoclave • Selecting/ starting the cycle via the touchscreen display • Cycling through the menus on the touchscreen display • Selecting drying and non drying modes via the touchscreen display • Cleaning the autoclave • Replacing filters and door seals • Downloading cycle data

Fig. 1 Main Components



- | | | | |
|---|--------------------------------------|----|---------------------------------|
| 1 | Chamber door | 8 | Serial Label |
| 2 | Door handle | 9 | Reservoir water drain |
| 3 | Door catch | 10 | Water reservoir lid/ steam vent |
| 4 | Touch screen display (Fig. 2) | 11 | Mains connection and fuses |
| 5 | Mains on/off switch | 12 | Door seal |
| 6 | RS232 port | 13 | Chamber |
| 7 | USB host port (for USB memory stick) | | |

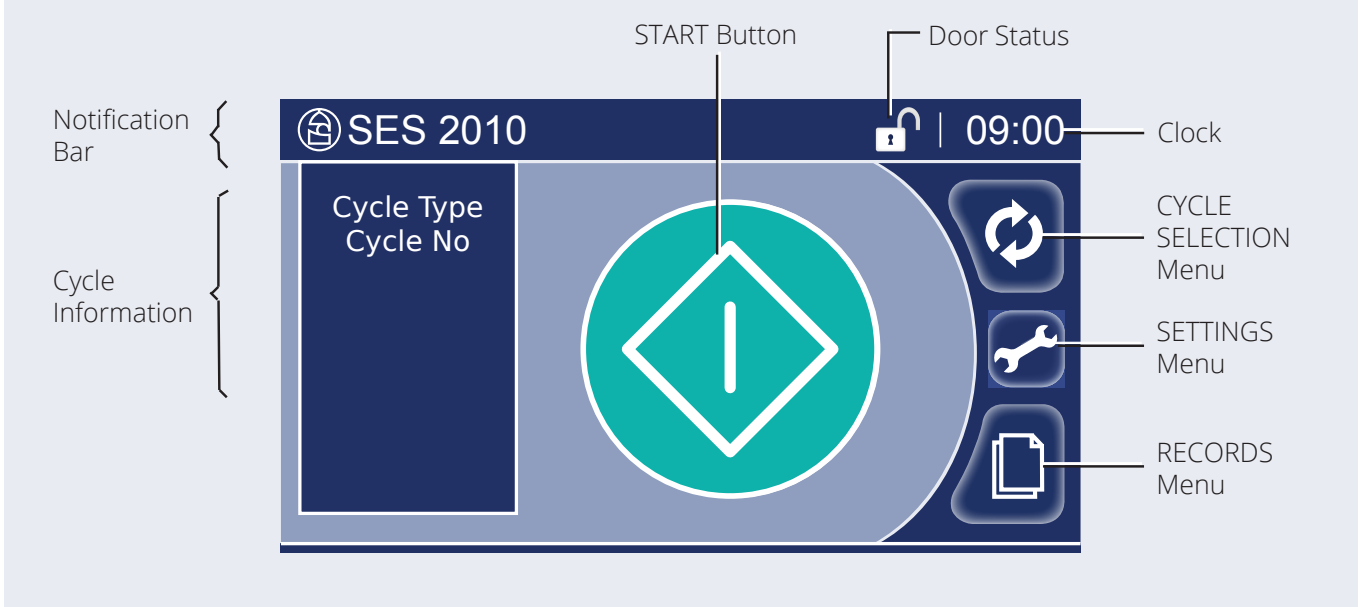Fig. 2 Typical START Screen



2 Introduction & installation

2.1 Introduction

The Eschmann SES 2010 is an 'N' type, non-vacuum autoclave with an 11 litre chamber. The autoclave runs a 134°C cycle with selectable drying.

The autoclave is intended for the sterilisation of solid, unwrapped instruments and has been designed in accordance with applicable standards.

2.2 Installation

We recommend that the autoclave is installed and commissioned by an Eschmann or Eschmann Trained Engineer. Contact Eschmann for installation, User Training and Warranty Registration.

Note: *This manual does not include configuration settings to be made only by an Eschmann or Eschmann Trained Engineer*

2.3 Unpacking

Remove the outer packaging and all packaging pieces. Recycle packaging in accordance with local recycling practices.

CAUTION

Always move or lift the autoclave using at least two people. Observe good manual handling techniques at all times, with hands placed under the base on each side of the autoclave.

Check for the following items (some are packed inside the chamber - remove these):

User Manual	Mains power cable
Drain tube	Tray carrier and tray lifter

Please note that contents may vary

2.4 Autoclave positioning

- Place the autoclave on a flat, clean, dry, level surface that is water resistant, heat resistant and suitable for the weight.
- Ensure that the mains plug and power switch will be easily accessible during use.
- There must be a clearance of at least 50mm to the sides and rear of the unit and free space above for filling the water reservoir.
- Clearance to the front to open the door: 300mm
- Ensure that the door is able to open fully.
- For autoclave environmental requirements, refer to Appendix 5.















2.5 Mains connection

Connect the autoclave to a 230V AC mains power socket using the supplied power cable.

Use only the Eschmann supplied mains power cable. It is fitted with the correct fuse.

Eschmann recommend the use of a suitable RCD (Residual Current Device).

2.6 Information/ safety symbols

Product Safety Symbols	
 or 	Attention: See the instructions for use
	Consult the instructions for use
	DO NOT use tap water
	Mains input fuse rating
	Potentially hot surfaces
Product Information Symbols	
	Manufacturer
	Date of manufacture
	Serial number
	Catalogue number
	WEEE recycling: Do not dispose of product at end of life with other waste
	The UKCA marking of the product certifies that it complies with the Medical Devices Regulations 2002 (SI 618), as amended by the EU Exit Regulations 2019 (SI 791) and 2020 (SI 1478)
	The CE marking of the product certifies that it complies with the Medical Devices Directive 93/42/EEC
	EC authorised representative

3 Controls

3.1 Touch screen & menu operation

Touch screen



The autoclave must be switched on to use the touch screen.


The autoclave uses a resistive touch screen. This requires a light pressure to operate. You can use a gloved or non-gloved hand.


Use the touch screen to view and navigate all the menus and pages contained in the menu system.


Menu operation

In the menus:

-  Use the up/down buttons to highlight your chosen setting in the blue area of the screen.
-  For example, pressing the up arrow will effectively move the blue area upwards.

 Press to accept a setting or selection

 Press to return without change

 Press to edit a selection

3.2 The Keypad



At points throughout the menu, whenever the keypad is displayed, you must enter the required information. For example; PIN code (Personal Identification Number), date, time, delay time, cycle number.

Each type of keypad expects the input to be in a certain format and will prevent incorrect input.

The formats are:

PIN code	1234	four digits
Date	dd/mm/yy	day/month/year
Time (24hr clock)	hh/mm	hours/minutes
Delay time	hh/mm	hours/minutes
Cycle number	123456	six digits

Below is the method for entering a four digit PIN code. The process is similar for all other keypads.

Type your four-digit PIN code:

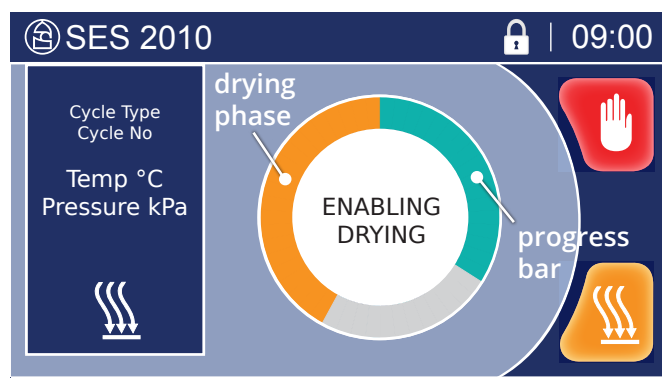
Each time a number is pressed, the flashing line (underscore) is replaced with the number typed. The line then highlights the next number.

- Press **OK** to accept the PIN code.
 - If the PIN code is accepted, the next screen will be displayed.
 - If the PIN code is not accepted the Autoclave beeps and starts the process of entering a PIN code again.
- Press Backspace (↵) to re-enter an incorrectly typed number. A final backspace will display the previous screen.

If your PIN code did not work:

- You may have entered it incorrectly.
- You may not have entered four digits. Pressing OK with less than four digits displayed will start the process of entering a PIN code again.
- The system may require a higher level of PIN code. For example, certain Menus and Warnings require a Manager PIN. However, the main process of running a cycle only ever requires a User PIN.

3.3 The Progress screen



The Progress Screen indicates cycle progress, cycle phase and cycle information (e.g. ambient temperature and pressure).

To do this, a green circular progress bar is drawn on the screen until it is completed for the selected cycle. At this point the cycle ends.

When drying is deselected for a cycle, the progress bar does not include the (orange) drying phase.

The phases of the cycle are indicated in the centre of the progress screen, as follows:

	Preconditioning Phase
	Sterilisation Phase
	Post Sterilisation Cooling Phase
	Drying Phase

Two-minute indication

The autoclave marks the two-minute point in the Sterilisation Phase for manual record keeping.

It displays the following, in order, to count down to the two minute mark:



The two-minute mark has been reached the moment **2** is displayed.

3.4 Notification screens

These may be displayed before and during a cycle.

Complete the actions requested to continue using the autoclave.

There are three types of notification screen:

Fill water reservoir	Section 4
Service screens	Section 10
Advisory screens	Section 11

3.5 Autoclave manager

Appoint a staff member to be the Autoclave Manager.

The Autoclave Manager, who must be trained in the use of the autoclave and will control autoclave use.

The Autoclave Manager is responsible for:

- Ensuring all users are regularly trained in the use and maintenance of the autoclave.
- Maintaining a record of training attendance including evidence of trainee understanding.

Refer to Section 13: "Security" to set up security PINS and choose a security level for the autoclave.

Note: *REMEMBER YOUR "MANAGER" PIN CODE! A forgotten Manager PIN code will require a Service call-out to access the menus.*

3.6 Water quality

The water reservoir must be filled with water of an appropriate quality that is low in dissolved solids and has a low microbial count - do not fill during a running cycle.

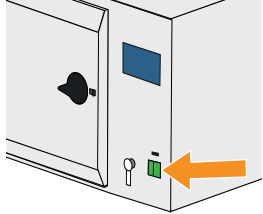

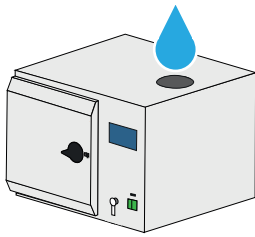


We recommend the use of:

- water treated by the reverse osmosis process
- distilled water
- deionised water

Refer to local authority requirements if in doubt.






- TDS (Total Dissolved Solids) to be less than 7mg/l (milligrams per litre), or 7ppm (parts per million)
- pH to be between 5.0 - 7.5
- conductivity <15 microSiemens/cm at 20°C

4 Preparing the autoclave

<p>1</p> 	<p>Switch the autoclave on</p>  <p>The Welcome screen is displayed. Note that the autoclave may take about one minute to initialise</p>	<p>Mains Switch : See Fig. 1, page 4.</p>
<p>2</p> 	<p>Fill the water reservoir</p>  <p>This screen is displayed if the water reservoir needs filling.</p> <ol style="list-style-type: none"> 1. Remove (Lift) water reservoir lid on the top of the autoclave. 2. Carefully pour in water to the 'MAX' mark. Refer to Section 3.6: "Water quality." 3. Replace reservoir lid by pushing into place. 	 <p>DO NOT use tap water!</p> <p>You can fill the water reservoir while the autoclave is switched on.</p> <p>DO NOT fill the water reservoir while running a cycle.</p> <p>The water reservoir MUST be emptied daily. Refer to Section 9.</p>

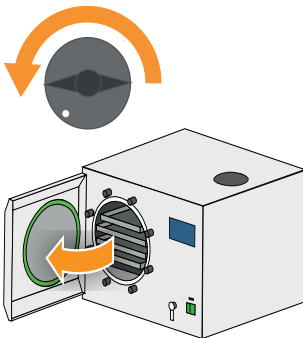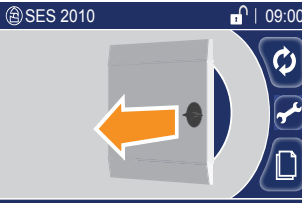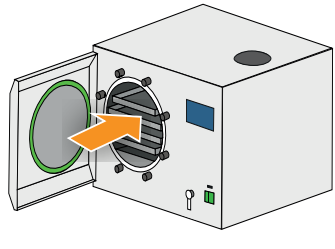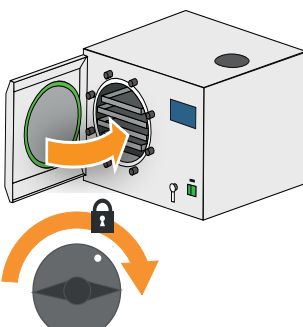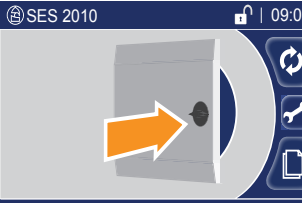⚠ WARNING!
 Be careful not to splash water on the top of the autoclave.
 If water is poured onto the autoclave, disconnect from the mains socket and dry thoroughly before use













5 Loading the autoclave

<p>1</p> 	<p>Open the door</p>  <ol style="list-style-type: none"> 1. Turn door handle anti-clockwise and release door catch. See Fig 1, pg 4. 2. Open the door. 	<p>Clean the door seal and chamber face every day using Eschmann recommended wipes.</p>				
<p>2</p> 	<p>Load the instruments to be sterilised</p> <table border="1" data-bbox="566 1512 1125 1668"> <thead> <tr> <th>Load</th> <th>Cycle type</th> </tr> </thead> <tbody> <tr> <td>Solid, unwrapped instruments intended for immediate use.</td> <td>N</td> </tr> </tbody> </table> <p>Maximum autoclave load: 4.5kg. Load trays using the supplied tray lifter.</p>	Load	Cycle type	Solid, unwrapped instruments intended for immediate use.	N	<p>All instruments MUST be thoroughly cleaned before sterilisation, ideally in an Eschmann Instrument Washer Disinfector.</p> <p>Refer to Appendix 2 for further loading information.</p>
Load	Cycle type					
Solid, unwrapped instruments intended for immediate use.	N					
<p>3</p> 	<p>Close the door</p>  <ol style="list-style-type: none"> 1. Push the door closed. 2. Turn door handle fully clockwise, so that the white dot is below the 'closed padlock' symbol. 	<p>Check there are no obstructions.</p> <p>Be careful not to trap your fingers when closing the door.</p>				

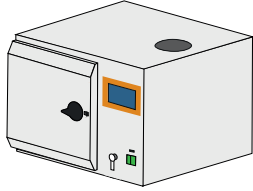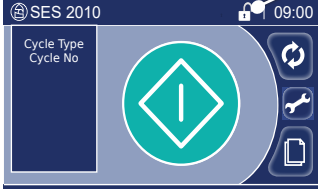6 Starting and running a cycle

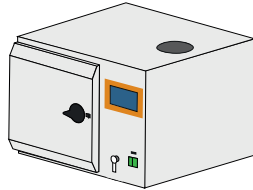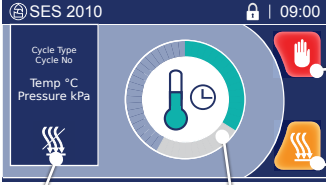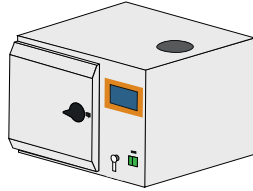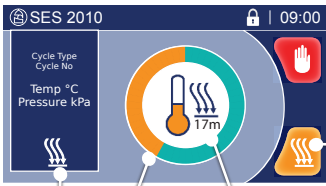Cycle Chart

Cycle type	Load type
134° N	Solid instruments (unwrapped) - with selectable drying.

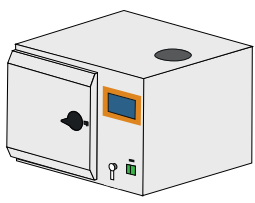
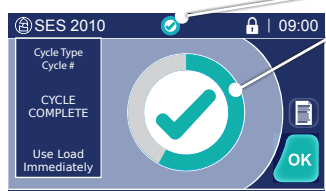

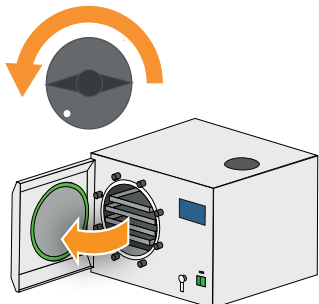
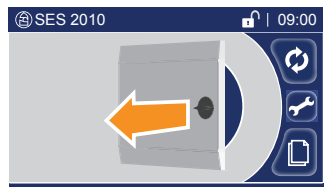
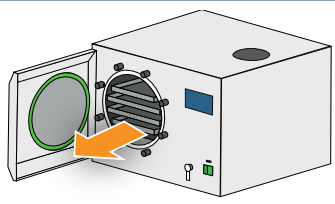
<p>1</p> 	<p>Start the cycle</p>  <p>Padlock on screen Indicates door lock status:  door locked  door unlocked</p> <p> 1. Press and hold for 1 second to start the cycle. 2. The door will now automatically lock and the cycle will start. The Cycle Progress Screen will be displayed.</p>	<p>You may be prompted to enter a PIN code. Refer to Section 3.2</p> <p> Note: If the autoclave too hot symbol is displayed, the autoclave will cool before starting the cycle.</p>
<p>2</p> 	<p>Cycle progress</p>  <p>Indicates drying not selected</p> <p>Cycle progress</p> <p>Cycle abort</p> <p>Drying selection button</p>	<p>Refer to Section 3.3 for description of the cycle progress screen</p> <p>Drying is off by default.</p>
<p>3</p> 	<p>Selecting/ de-selecting drying</p>  <p>Press the drying selection button during the cycle to turn drying on or off.</p> <p>Drying progress</p> <p>Drying time remaining</p> <p> Indicates drying selected</p> <p> Indicates drying not selected</p>	<p>The drying portion of the progress bar is orange when drying is selected, and grey when drying is not selected.</p> <p>Drying time: 17 minutes</p>

WARNING!







During operation:
 Do not remove or cover the water reservoir lid (Fig 1, pg 4).
 Caution: hot steam is vented.

7 Cycle completion and unloading


WITHOUT DRYING

<p>1</p> 	<p>Cycle complete</p>  <p>Green tick indicates cycle passed.</p> <p>OK Press to unlock door when OK button displayed.</p>	<p>You will be prompted to enter a PIN code. Refer to Section 3.2.</p>  Press to view the completed cycle log
<p>2</p> 	<p>Open the door</p>  <ol style="list-style-type: none"> 1. Turn door handle anti-clockwise and release door catch. See Fig 1, pg 4. 2. Open the door. 	<p>If an action is required, for example, the Fill Water Reservoir screen is displayed, you can still open the door.</p>
<p>3</p> 	<p>Unloading</p> <p>Use the tray lifter - refer to Appendix 2.</p> <p>THE UNIT IS NOW READY TO RUN ANOTHER CYCLE.</p> <p>Leave the door ajar between cycles.</p>	<p>Take care! The load, chamber and door will be VERY HOT.</p> <p>Sterilised items are intended for immediate use.</p>

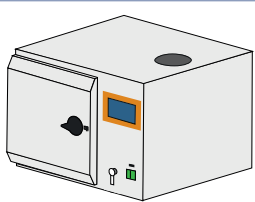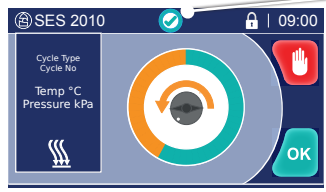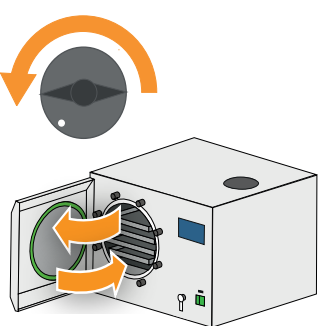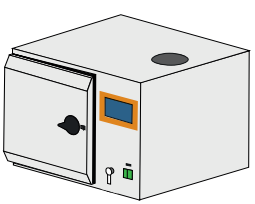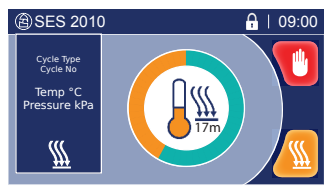
WITH DRYING

<p>1</p> 	<p>Sterilisation complete</p>  <p>Green tick indicates cycle passed.</p> <p>OK Press OK to unlock door, when OK displayed. You will be prompted to enter a PIN code. Refer to Section 3.2.</p>	<p>Green tick indicates cycle passed.</p>
<p>2</p> 	<p>Positioning door for drying</p> <ol style="list-style-type: none"> 1. Turn door handle anti-clockwise and release door catch. See Fig 1, page 4. 2. Open the door briefly to allow residual moisture to escape. 3. Push the door closed again so that it clicks onto the catch, slightly ajar. 	<p>Take care! The load, chamber and door will be VERY HOT.</p>
<p>3</p> 	<p>Drying phase</p>  <p>Drying will run for 17 minutes.</p>  Press to end drying early. <p>Door on catch, slightly ajar.</p>	

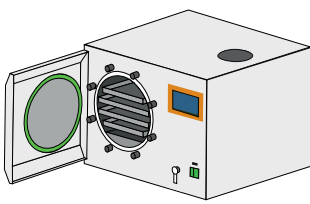4



Drying Completed

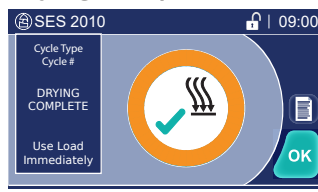


OK Press to accept cycle.
You will be prompted to enter a PIN code. Refer to Section 3.2.


Open the door by releasing the catch. Unload autoclave following Section 7, step 3.

Press to view the completed cycle log.

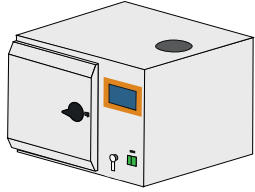5



If door NOT opened for drying

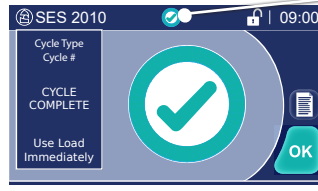

OK Press to accept cycle.
You will be prompted to enter a PIN code. Refer to Section 3.2.

Unload autoclave following Section 7, steps 2 and 3

Green tick indicates cycle passed

This screen is only displayed if drying was selected, but the door was not opened onto the catch.

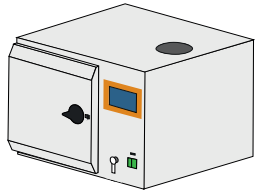
8 How to stop the cycle

⚠ WARNING!


 Press and hold for 2 seconds to stop the cycle.


The load is NOT sterilised.
When a sterilising cycle is aborted the load must be treated as non-sterile. The load must be sterilised by running the cycle again.

1




Stop the cycle

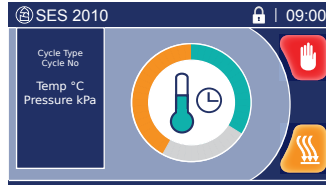


 Press and hold for 2 seconds to stop the cycle.

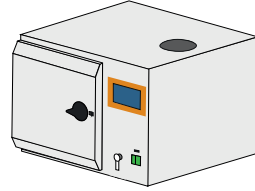2



The cycle has stopped

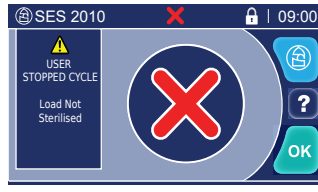


Red cross indicates cycle failed.

OK

When OK appears after a short delay, you will be prompted to enter a PIN code. Refer to Section 3.2.

There may be a delay whilst the chamber stabilises, allowing the door to be opened.

 Press this to display details of the cycle.
Refer to Section 11 for further details.

Take care!
The load, chamber and door will be VERY HOT.

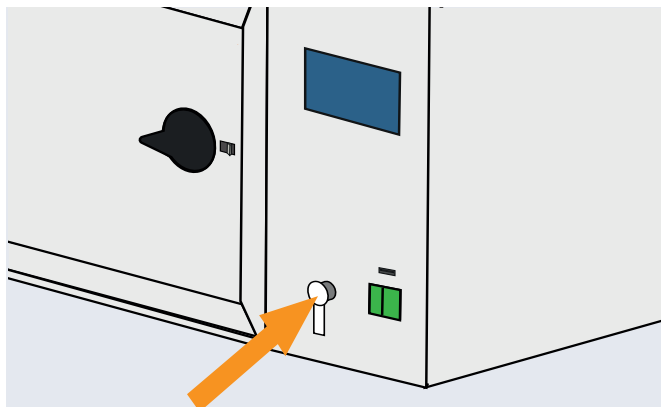
9 After use

9.1 Emptying the water reservoir

We recommend that the water reservoir is emptied, as a minimum, at the end of each day.

⚠ WARNING!
Water in the reservoir can be very hot, take care when draining to avoid contact with the hot water.

The water reservoir drain point is located as shown below:




Remove the drain plug from the water reservoir connector by pressing the metal thumb release on the connector side.

Place the cut end of the drain tube into the sink, or a container on a surface below the level of the autoclave (Capacity: 5 litres or more).

Connect the drain tube to the water connector by pressing the metal thumb release on the connector side. This action opens the valve. The red drain tube clip is now holding the water back.

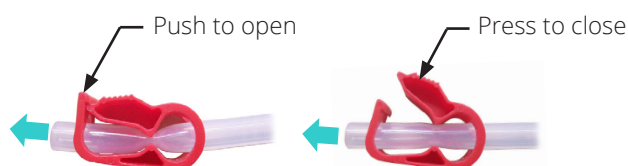Open the drain tube clip to drain the water.



If using a container, remember to empty it. A five litre container is sufficient to empty reservoir.

Once drained, release the drain tube from the connector by pressing the thumb release on the side of the connector. Refit the drain plug.

Check that the drain tube clip is not damaged or leaking. If so, **replace as soon as possible** - refer to Appendix 5: "Technical data".

9.2 Switch off

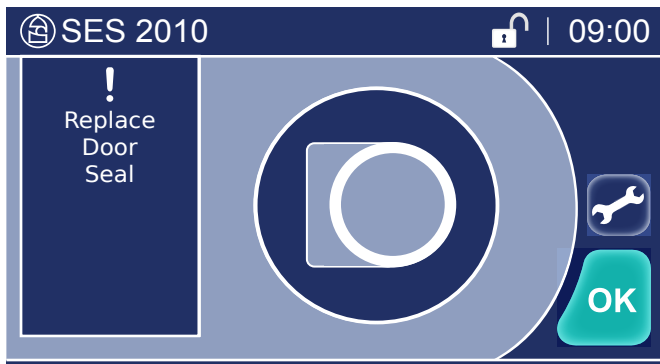
Switch off the autoclave at the mains switch (Fig 1, pg 4) when not in use. Leave the door ajar and the chamber empty.

9.3 Cleaning

At the end of each day carry out the Cleaning and Care instructions. Refer to Appendix 3.

10 Service screens

10.1 Replace the door seal



The door seal needs replacing when this screen is displayed.

You can replace the door seal with the autoclave switched on.

OK Press OK to temporarily clear the screen.



WARNING!

The water, chamber and load could be very hot and potentially non-sterile.

1. The door seal should be replaced by a qualified Eschmann Engineer. Alternatively a user may replace the door seal through use of the Door Seal Kit, available from Eschmann. The kit comprises comprehensive instructions and all required parts to ensure safe and correct replacement. Part No. MUW216747.
2. Follow local biohazard procedures for cleaning and decontamination when disposing of the door seal.
3. Update the Settings Menu counter for the door seal.

Updating the Settings Menu

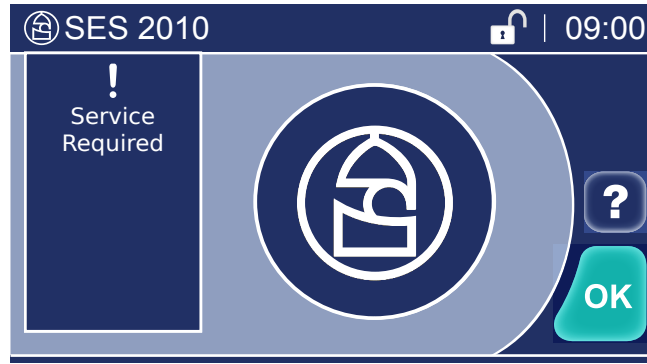
It is important to update the door seal counter in the Maintenance Menu. The counter counts the number of cycles performed using the seal. This allows the autoclave to warn you when replacement is required.



Press to display the Settings Menu from the Open Door, Close Door, or Cycle Start screen and select Maintenance.

Refer to Section 12.2 : 'Settings Menu' - (Managers only).

10.2 Service required



The autoclave requires servicing by a qualified Eschmann, or Eschmann trained engineer. Contact Eschmann as soon as possible.



Press to display the details required for contacting Eschmann.



Press to acknowledge the warning and clear the screen. The screen will reappear at the beginning of each cycle until the autoclave is serviced.

11 Advisory screens

Advisory screens may be displayed at any time to protect the autoclave.

When an Advisory screen is first displayed some touch screen buttons are not available whilst the autoclave is depressurising in preparation for the door being opened.

This may take up to two minutes, at which point one of the following advisory screens will be displayed.

Refer to Appendix 4. Trouble shooting, Table 1.

CAUTION

The water, chamber and load could be very hot and potentially non-sterile.

After clearing any warning, always allow the chamber to dry before starting another cycle.

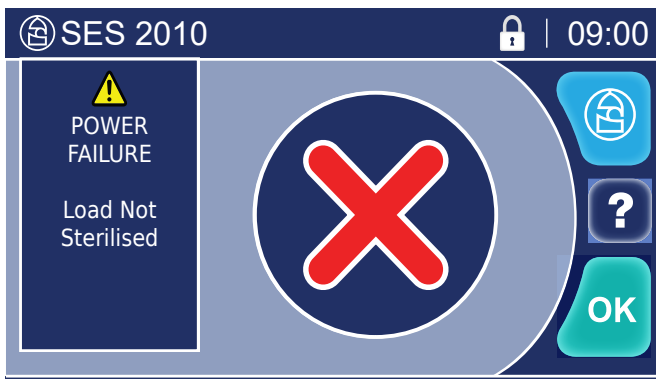
⚠ WARNING!

If an advisory screen is displayed, the load must be treated as **NON-STERILE**.

Sterilise the load by running the cycle again.

Do not attempt to solve the problem using tools or by tampering with the autoclave.

POWER FAILURE



Displayed if the autoclave experiences a power interruption. Note: Door lock will stay in (Or change to) the locked state, upon power failure.



Press to display the details required for contacting Eschmann.

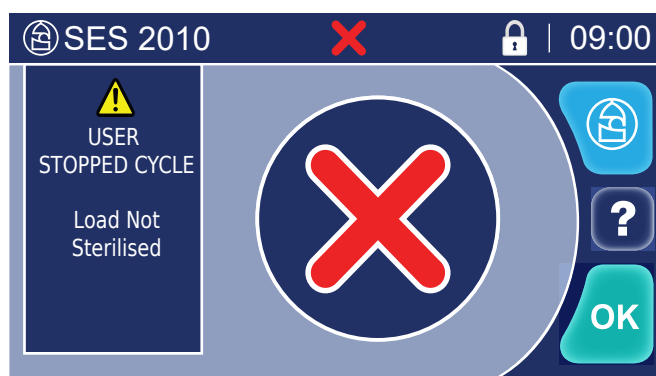


1. When OK appears after a short delay, press to acknowledge the warning and clear the screen.
2. Check mains supply, mains connections, autoclave fuses and mains cable fuse.
3. Run the cycle again. If the warning persists, contact Eschmann.



Eschmann use only (PIN protected).

USER STOPPED CYCLE



Displayed when the STOP button  is pressed to stop a cycle.



Press to display the details required for contacting Eschmann.

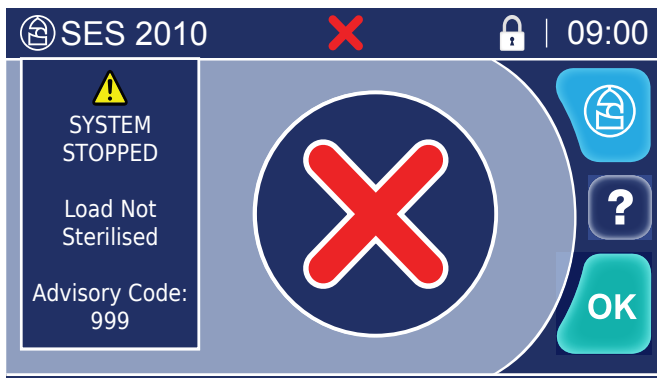


When OK appears after a short delay, press to acknowledge the warning and clear the screen.



Eschmann use only (PIN protected).

SYSTEM STOPPED



Displayed when the autoclave prevents operation **DURING A CYCLE**.



Press to display the details required for contacting Eschmann.

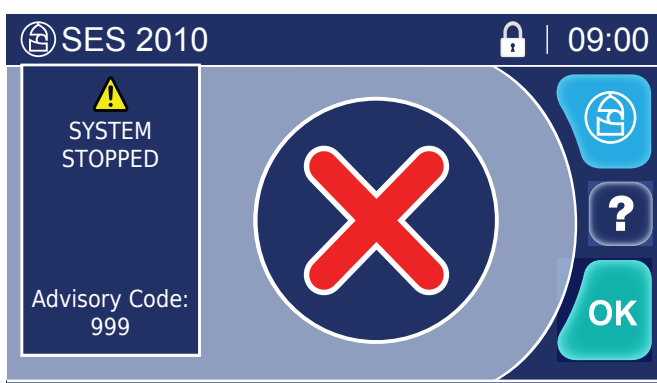


1. Note the Advisory Code. When OK appears after a short delay, press to acknowledge the warning and clear the screen.
2. Refer to Appendix 4 and complete the Actions against the relevant Advisory Code.
3. Run the cycle again. If the warning persists, contact Eschmann.



Eschmann use only (PIN protected).

SYSTEM STOPPED



Displayed when the autoclave prevents operation when a **CYCLE IS NOT IN PROGRESS**



Press to display the details required for contacting Eschmann.



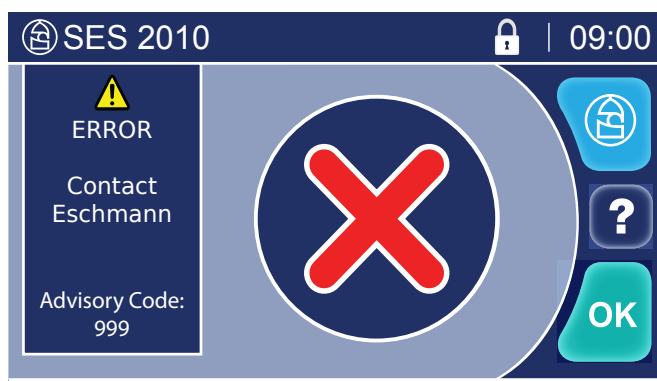
1. Note the Advisory Code. When OK appears after a short delay, press to acknowledge the warning and clear the screen.
2. Refer to Appendix 4 and complete the Actions against the relevant Advisory Code.
3. Run the cycle again. If the warning persists, contact Eschmann.



Eschmann use only (PIN protected).

ERROR

All errors require you to contact Eschmann.




Displayed when the autoclave experiences an Error.



1. When OK appears after a short delay, press to acknowledge the warning and unlock the door.
2. The door locked symbol at the top of the screen will change to the door unlocked symbol when the door is ready to be opened to retrieve instruments.

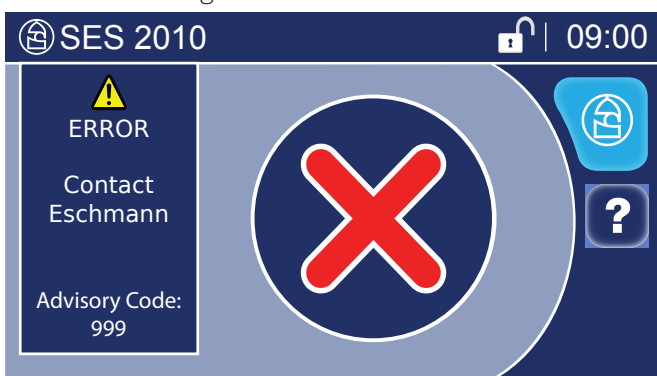The screen changes to:



The autoclave is prevented from being used as this Error screen cannot be cleared.



Press to display the details required for contacting Eschmann.

Switch the autoclave off and unplug from the electrical wall socket.

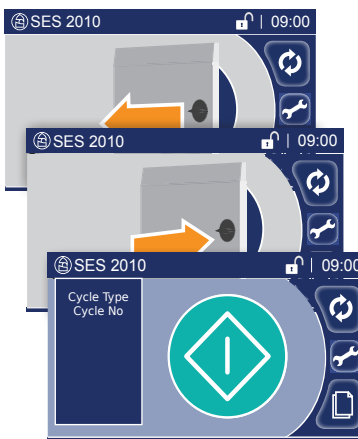

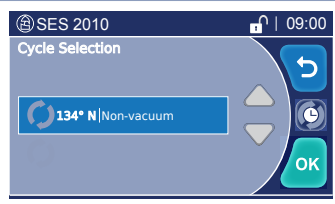






Eschmann use only (PIN protected).

12 Menus

12.1 Cycle selection menu

Selecting a cycle







1		<p>The cycle has not yet started</p> <p>From the Open Door, Close Door, or Cycle Start screens:</p>  Press to view the Cycle Selection menu.	<p>The 134° N cycle is the factory default for each new cycle.</p> <p>The current cycle type is shown on the Start screen.</p>
2		<p>Select the cycle</p>  Use the up/down buttons to highlight your chosen cycle in the blue area of the screen.   Press to accept and return to your starting screen in step 1 above.	 Press to return without change.

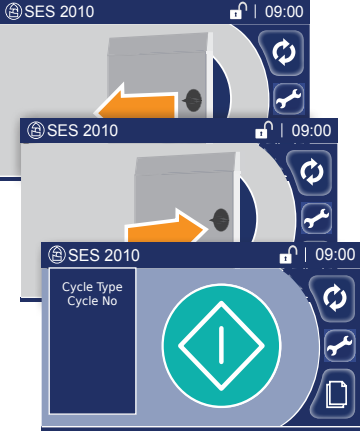Setting a delayed start

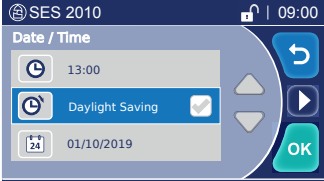





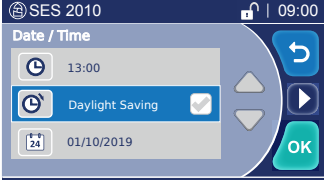


<p>1</p> 	<p>The cycle has not yet started The Autoclave is loaded</p> <p>From the Open Door, Close Door, or Cycle Start screens:</p> <p> Press to view the Cycle Selection menu.</p> <p>Delayed Start allows you to delay the start of a cycle by up to 23 hours and 59 minutes.</p>	
<p>2</p> 	<p>Select the cycle to delay</p> <p> Use the up/down buttons to highlight your chosen cycle in the blue area of the screen.</p> <p></p> <p> Do not press.</p> <p> Press to view the Delay Start screen.</p>	<p> Press to return without change.</p>
<p>3</p> 	<p>Set a delayed start time</p> <p> Press to accept the delay time shown and go to 5 below or</p> <p> Press to enter a new delay time.</p>	<p> Press to return without change.</p>
<p>4</p> 	<p>Insert the Delay Time</p> <p>Enter the delay time in hours and minutes using the 24 hour clock.</p> <p>Press OK to accept.</p>	<p>The maximum delay time is 23 hours and 59 minutes.</p> <p> Press to return without change.</p>
<p>5</p> 	<p>Start the Delayed Cycle</p> <p> Press to start the delay cycle. The door will lock and the delayed cycle will start.</p>	<p>The delay time begins when the Start button is pressed.</p>
<p>6</p> 	<p>Delayed Cycle Running</p> <p>When the count down is complete (0:00), the Run Screen will be displayed for the duration of the cycle.</p> <p> Press to stop the cycle.</p> <p> Toggle to select drying ON/OFF</p>	<p>A count down is shown on the screen.</p> <p>If drying is selected, the door will need to be opened onto the catch upon completion of the sterilisation phase. Refer to Section 7.</p>

12.2 Settings Menu (Managers only)


Viewing the Settings Menu

<p>1</p> 	<p>The cycle has not yet started</p> <p>From the Open Door, Close Door, or Cycle Start screens:</p>  Press to view the Settings menu.	<p>The Settings menu contains:</p> <ul style="list-style-type: none">Date / TimeSoundsPIN ManagementSecurity LevelServiceMaintenanceAbout SES 2010
<p>2</p> 	<p>Select the required setting</p>  Use the up/down buttons to highlight your chosen setting in the blue area of the screen.  Press to view the setting screen.	 Press to return without change.


Date / time

<p>1</p> 	<p>Set the autoclave's date and time</p>  Use the up/down buttons to highlight your chosen setting in the blue area of the screen.  Press to edit the highlighted setting: Date or Time.	 Press to return without change.	<p>You must set the correct date and time as this is logged by the cycle records.</p>
<p>2</p> 	<p>Insert the date or time</p> <p>Enter date or time in the relevant keypad. Press OK when complete.</p> <p>The date format is Day/Month/Year. The time setting uses the 24 hour clock.</p>	 Press to return without change.	
<p>3</p> 	<p>Set the autoclave's date and time</p>  Press OK to accept.	 Press to return without change.	

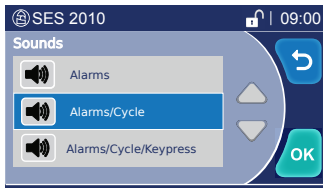Sounds



Select the sounds

 Use the up/down buttons to highlight your chosen setting in blue.



 Press **OK** to accept.

Alarms	Warning screens are announced by 5 short beeps
Cycle	Cycle start: short beep Cycle end: long beep
Keypress	Short beep with each key press



Press to return without change.

There are three sound types to choose from:

Alarms
Alarms/Cycle
Alarms/Cycle/Keypress

Note that alarms are always sounded.

The selection highlighted in blue is the current selection.

PIN management



Setting a PIN Code

Refer to Section 13.3: 'Setting a PIN code'.



Press to return without change.

Security level



Selecting the Security Level

Refer to Section 13.4: 'Selecting the security level'.



Press to return without change.

Service



Service

The Service menu is protected by a Service PIN Code that is for Eschmann use only.



Press to return without change.

Maintenance



Resetting the Door Seal Counters

The autoclave counts the number of cycles performed by the Door Seal. From this, it determines when the item will need to be renewed. The autoclave uses a simple traffic light system to indicate the usage of the parts:



Possible cycle quality impact - replacement overdue.



Immediate replacement required.



Cycle count not high enough to require replacement, although if the part is damaged, replace immediately.



Use the up/down buttons to highlight the chosen selection in blue.



Press to reset the highlighted selection.



Press **OK** to accept.



Press to return without change.

Replace the door seal before resetting the counter.



When you reset an item on this screen, it is marked with an arrow to indicate that you have selected to reset it.

Item(s) are reset when you press **OK**.

The **Replace the Door Seal** service screen begins to display when an amber warning is in force.

Configure Cycles

<p>1</p> 	<h3>Enable Cycles</h3> <p>Use this screen to enable the cycles available to the user on the Cycle Selection screen.</p> <ul style="list-style-type: none"> Use the up/down buttons to highlight the chosen selection in blue. Press to toggle the check box. A green tick indicates the cycle type is enabled. A greyed out tick indicates the cycle is not enabled. Press OK to accept. The default cycles screen (See below) will be displayed.	 Press to return without change. Currently the SES 2010 has a 134°C N type cycle only.
<p>2</p> 	<h3>Default Cycles</h3> <p>Use this screen to select the autoclave default cycle. Only one default cycle can be selected.</p> <ul style="list-style-type: none"> Use the up/down buttons to highlight the chosen selection in blue. Press to toggle the check box. A green tick indicates the cycle is selected as the default cycle. A greyed out tick indicates the cycle is not selected as the default cycle. Press OK to accept.	 Press to return without change.
<h2>About SES 2010</h2>		
	<h3>About SES 2010</h3> <p>Contains useful information about your autoclave when contacting Eschmann.</p> <ul style="list-style-type: none"> Press to return to the starting screen. Press to view the Records menu.	The last advisory code may also be presented on this screen.

12.3 Cycle Records

Cycle records can be downloaded onto a USB memory stick. Follow the 'Downloading Records' instructions (1 to 5) in this section. (Note: The USB memory stick should be <16Gbs USB FAT32 format only).


If a USB memory stick is left in the USB port (Fig. 1, item 7), the autoclave will download each cycle record automatically upon cycle completion.

Cycle records can also be viewed directly on the screen. Follow the 'View Cycle' instructions in this section.

Downloading Records

<p>1</p> 	<p>From the Open Door, Close Door, or Cycle Start screens:</p> <p> Press to view the Records menu.</p> <p>Insert a USB memory stick into the USB port. (Fig. 1, item 7)</p>	
<p>2</p> 	<p>Download records</p> <p> Use the up/down buttons to highlight the chosen selection in blue.</p> <p></p> <p> Press to accept and view the next screen.</p>	<p>Select from:</p> <ul style="list-style-type: none"> Successful Cycles All Cycles Failed Cycles <p> Press to return without change.</p>
<p>3</p> 	<p>Download records</p> <p> Use the up/down buttons to highlight the chosen selection in blue.</p> <p></p> <p> Press to accept and view the next screen.</p> <p>All Since Last Download will start download of all cycle records since the last download.</p> <p>See below when selecting By Cycle Number and By Date.</p>	<p>Select from:</p> <ul style="list-style-type: none"> By Date All since last Download By Cycle Number. <p> Press to return without change.</p>
<p>- By Cycle Number</p>		
	<p>Set FROM / TO Cycle Numbers</p> <p> Use the up/down buttons to highlight the chosen selection in blue.</p> <p></p> <p> Press to enter the selected cycle number (See below).</p> <p> Press to accept the cycle numbers shown and begin download.</p>	<p> Press to return without change.</p> <p>Input boxes default to the most recent cycle number, but can be changed.</p>
	<p>Enter the cycle number</p> <p>Enter the cycle number in the keypad. Zero and numbers prefixed by zeros are not accepted.</p> <p>Press OK to accept the cycle number entered.</p>	<p> Press to return without change.</p>

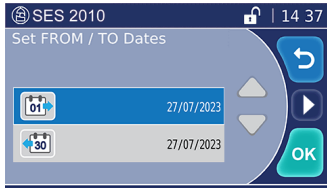- By Date



Set FROM / TO Dates

Use the up/down buttons to highlight the chosen selection in blue.



Press to enter the selected date. (See below)

Press to accept the dates shown and begin download.



Press to return without change.

Input boxes default to today's date, but can be changed.



Enter the date

Enter the date in the keypad.

Press **OK** to accept the date entered.

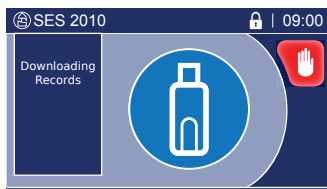
The date format is Day/Month/Year.



Press to return without change.

Download

4



Downloading Records

The records are being downloaded via the USB port.

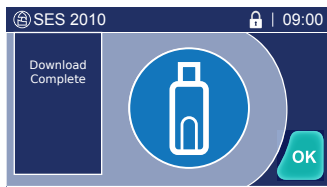


Press to stop the download before download is complete.

If **Plug In USB** is displayed, check that a USB memory stick is correctly connected.

Do NOT remove the USB memory stick during download.

5



Downloading Complete

Download is complete. Exit the Records menu before removing the USB device.



Press **OK**

If **Download Failed** is displayed, check for the following:

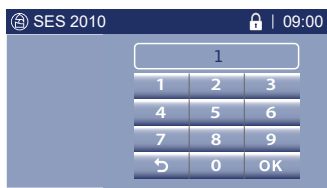
Memory stick not formatted

Memory stick is full

Memory stick not plugged in

Search returned no data

View Cycle



To view a cycle record directly on the screen:

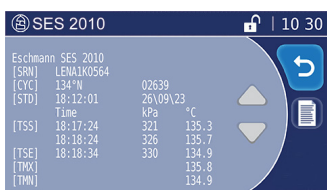
From the Download Records menu (See instructions 1 and 2 in this section), select 'View Cycle'.

Enter the cycle number in the keypad. Zero and numbers prefixed by zeros are not accepted.

Press **OK** to accept the cycle number entered.



Press to return without change.



The cycle record can now be viewed on the screen.



Press to view the detailed cycle record



Press If a printer is enabled, to obtain a printout of the cycle record.

13 User security

13.1 Overview

The autoclave can be configured to provide varying levels of security.

- This takes the form of a four digit PIN code which must be entered into the Keypad at different points in the touch screen menu.



Levels of security

The autoclave provides for three types of user:

Type	Description
User	Day-to-day user
Manager	The Autoclave Manager (more than one if required)
Service	Eschmann use only

- A User PIN will only work at a User location.
- A Manager PIN will work at both a User or Manager location.

There are fixed points in the touch screen menu that require you to enter a PIN code that cannot be disabled:

Fixed PIN Locations	Type
End of cycle	User/Manager
All Warnings	User/Manager
Service Required screen	User/Manager
Settings Menu	Manager
Records Menu	Manager

Default setting

The autoclave has twelve PINs. PIN 0001 is configured as a Manager allowing full access to the autoclave (except Eschmann Service screens). For security, this number should be changed by the Manager when first receiving the unit.

Note: REMEMBER YOUR "MANAGER" PIN CODE! A forgotten Manager PIN code will require a Service call-out to access the menus.

All other PINs are disabled by being set to 0000.

Security PINs are configured in the Settings\PIN Management menu. The Settings menu is protected by the Manager PIN.

13.2 Autoclave manager

Appoint a staff member to be the Autoclave Manager.

The Autoclave Manager, who must be trained in the use of the autoclave, will control the use of the Autoclave.

The Autoclave Manager is responsible for ensuring that the correct sterilisation cycle is used for the items being sterilised, ensuring that records are kept correctly, and for keeping the PIN code(s) confidential.

We recommend that you set up the security so that you have:

- 1 Autoclave Manager (perhaps more in large departments).
- 1 or more Users.

Reasons for security

BS EN 13060 specifies the general requirements and test methods for small steam sterilisers and sterilisation cycles used for medical purposes, or for materials that are likely to come into contact with blood or body fluids.

The choice of sterilisation cycle must be appropriate for a particular load type.

HTM guidelines state that:

- Permanent records of every sterilisation cycle should be kept.
- Written records of all testing and maintenance should be kept
- Sterilisation performance must be checked frequently.

The security offered by the SES 2010 Autoclave captures all of these elements. The PIN codes that you create form a major part of the record keeping for this Autoclave. They appear in the cycle records (Section 12.3 : "Records Menu"), and in records printed using the Cycle Logger application (Section 14).

The PIN codes will allow you to know:

Who started a cycle	High security level
Who stopped a cycle	All security levels
Who unlocked the door at the end of a cycle	All security levels
Who cleared Service screens	Medium and High security levels
Who needs training	All security levels

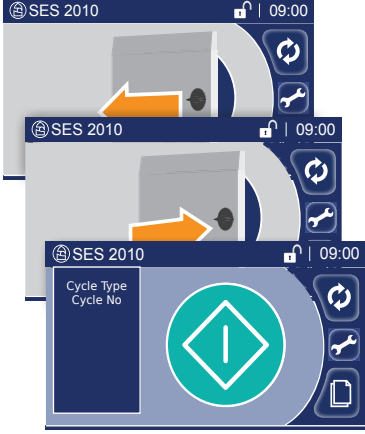

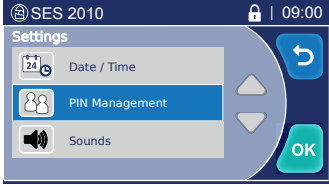

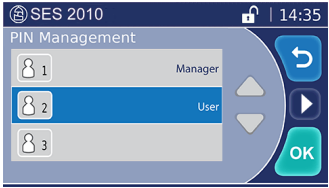






13.3 Setting a PIN code

To set a PIN	Enter a four digit number (except 0000)
To disable a PIN	Enter 0000

The PIN won't set?


If you type in a PIN but the screen reverts to 0000, then that PIN is already in use. Enter a different four digit number.

Note: *Eschmann reserve a number of PINS.*

<p>1</p> 	<p>The cycle has not yet started</p> <p>From the Open Door, Close Door, or Cycle Start screens:</p> <ul style="list-style-type: none">  Press to view the Settings menu. Manager PIN required to proceed. 	
<p>2</p> 	<p>Select PIN Management</p> <ul style="list-style-type: none">  Use the up/down buttons to highlight PIN Management in the blue area of the screen.  Press to accept. 	<ul style="list-style-type: none">  Press to return without change.
<p>3</p> 	<p>Select the PIN to edit</p> <ul style="list-style-type: none">  Use the up/down buttons to highlight your chosen PIN in blue area of the screen.  Press to edit PIN or status level. 	<ul style="list-style-type: none">  Press to return without change.
<p>4</p> 	<p>Select to Edit PIN or status level</p> <ul style="list-style-type: none">  Use the up/down buttons to highlight the PIN or status level in blue.  Press to edit: toggles between User/ Manager, or displays a PIN code keypad. <p>Enter a new PIN code to replace the PIN displayed and press OK to accept.</p>	<ul style="list-style-type: none">  Press to return without change.  Press to return without change.
<p>5</p> 	<p>PIN management</p> <ul style="list-style-type: none">  Press to accept and save all changes. 	<ul style="list-style-type: none">  Press to return without change.



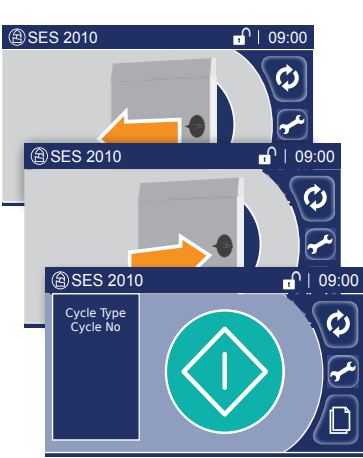13.4 Selecting the security level


- 1**






The cycle has not yet started

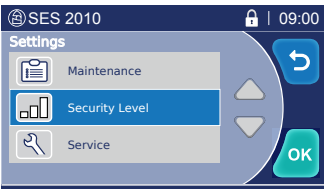
From the Open Door, Close Door, or Cycle Start screens:


 -  Press to view the Settings menu.
 -  Manager PIN required to proceed.
- 2**






Select Security Level


 -  Use the up/down buttons to highlight Security Level in the blue area of the screen.
 - 
 -  Press to accept.

 Press to return without change.
- 3**



Select required security level

 -  Use the up/down buttons to highlight the required Security Level in the blue area of the screen.
 - 
 -  Press to accept and return to the starting screen in 13.4, step 1 above.

 Press to return without change. MEDIUM security is recommended.

USER PIN locations enabled by selecting LOW, MEDIUM and HIGH Security Levels					
Security level	Service screens	Cycle select	Cycle start	Drying	Cycle end
LOW	✓				✓
MEDIUM Recommended	✓				✓
HIGH	✓	✓	✓		✓

14 Cycle Logger application

The digital storage of cycle data replaces the need for a traditional paper printer.

To view and print the files, transfer the files on to a PC using the USB memory stick.

(Note: The USB memory stick should be <16Gbs USB FAT32 format only).

14.1 Install the software

Suitable for use with Windows XP, 2000, Vista, and 7. Requires one USB port.

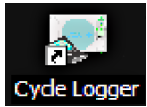
Disconnect the PC from the internet and any local area network (to prevent Windows® from installing the incorrect drivers).

Ensure there is no USB connection between the PC and the autoclave.

1. Insert the supplied USB memory stick that contains the Cycle Logger Application Software in to your PC. Follow the Installation Wizard's instructions.

If an older version is found, it will be uninstalled first. After this, remove and replace the USB memory stick to restart the installation.

2. When the Wizard has finished, remove the USB memory stick and shut down your PC.
3. Restart your PC.




The installation process is complete.

There is now a shortcut on your desktop to open the software.

14.2 Create a PC back-up

On a regular basis, (e.g. weekly), it is sensible to back-up your data from the USB memory stick to your PC, putting it into a back-up folder on your PC desktop.

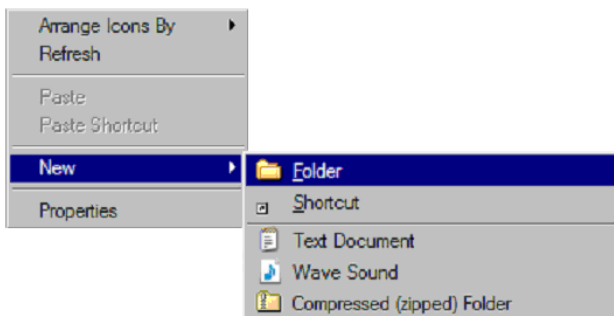
To create the back-up folder, right-click on the desktop and select *New*, then *Folder*.




Type in the name of the new folder. We suggest **SES 2010**. If you are backing-up more than one unit, you may want to add the unit's Serial Number to the title as well.

To back-up the files:

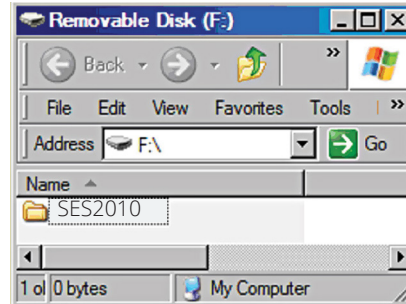1. Insert the USB memory stick directly into the USB port on your PC.
2. A window will appear saying AUTOCLAVE.
3. Double-click the AUTOCLAVE folder to display a folder similar to the one shown below. It will be named with the unit's Serial Number.



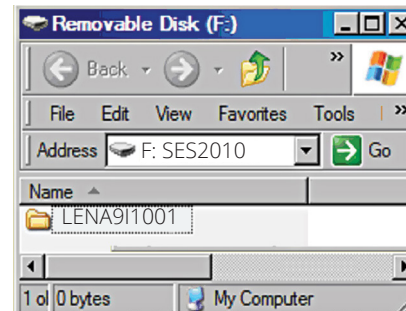4. To back-up the files, simply drag the folder onto your back-up folder on the desktop.



You are copying files in the same way you might copy photos from your camera. Partially minimise the windows on the PC to make drag-and-drop easier to use.

14.3 Viewing files

To view files in your back-up folder, use Windows® Explorer to view the folder contents.

Open single or multiple files by double-clicking any **.ELG** file.

Each file opens in its own window. Cycle files are ".ELG" files. They are named using the cycle number (e.g. 00186.ELG), or a code that is generated from the date and time (e.g. 1435874.ELG).

*If Windows does not recognise the file type, you can right-click the file and select *Open With*. Select the Cycle Logger program from the list. It will then recognise all **.ELG** files in the future.*

Viewing options

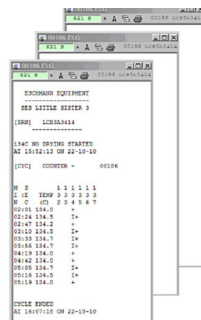


- A click to reduce the text size
- A click to enlarge the text size
- toggle to display all the data in the cycle, or only major events

Application Software

When an .ELG file is viewed, the Application Software is launched. It reads .ELG files and provides the following shortcuts, but you don't need to use it other than that.

Use these shortcuts when you have more than one window open:



List the open windows by right clicking in the Application Window. The open windows are listed in the Logs In Memory screen.

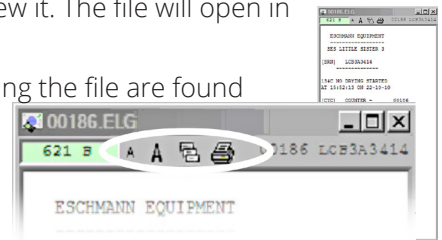
Bring a window to the front in the Logs In Memory screen by double-clicking its filename. Double-click it again to move the window to the top left of the PC screen.




14.4 Printing files

Select a file and view it. The file will open in its own window.

The tools for printing the file are found at the top of this window, as shown.



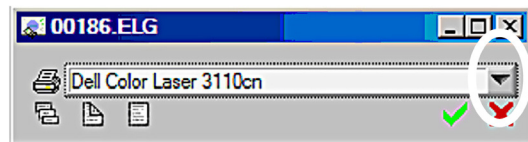1. Click to display the print options screen below.
2. Click the arrow to select your printer.



3. Click to print the file.

Printing options



- toggle to display all the data in the cycle, or only major events
- toggle to select double-sided or single-sided printing
- toggle to select portrait or landscape format
- print the file
- cancel printing

Fault finding

Application Not Loading

You must have Administrator rights to load programs. This is accessed on your PC via Start/Control Panel/User Accounts.

Compatibility Mode

To change the compatibility mode of the software, select Start/All Programs/Eschmann Equipment/Cycle Logger, and right-click on the second Cycle Logger.

Click on Properties. In the window click on the Compatibility tab. Check the Compatibility Mode check-box.

Select Windows XP from the drop-down menu. Click Apply and OK.

If problems still persist and you contact Eschmann we will need to know about your computer system.

We will need to know the:

- Windows version*
- Service Pack*
- Processor type*

This is identified by clicking the Start button on your PC. Right-click on My Computer (or Computer). Select Properties. Select the General tab to view the information.

Cycle output

The cycle output (opposite) can provide the following information:

[SRN]	autoclave serial number
[CYC]	cycle type (e.g 134°N) and cycle count (e.g. 00022)
[STD]	start time and date (hh:mm:ss, dd\mm\yy)

Sterilisation information (time, pressure and temperature every 60 seconds during sterilisation)

[TPP]	time of pressure point
[TSS]	time of sterilisation start
[p]	print sterilisation line ("all data" viewing only)
[TSE]	time of sterilisation end
[TMX]	maximum sterilisation temp
[TMN]	minimum sterilisation temp
[TDS]	time drying start
[TDE]	time drying end
[PTD]	time cycle finished processing ("all data" viewing only) (hh:mm:ss, dd\mm\yy)
[ETD]	end time and date (hh:mm:ss, dd\mm\yy)
[ERR]	error
[STS]	status at cycle end
[PCV]	pressure change value

Service and certification tags

[SVC]	service expired (date/cycle count of service expiry)
[CRT]	certification expired (date of certification)

Security log additions

(displays the number of the PIN: 1-12, 0 = not entered)

[DCB]	drying cancelled by
[CSB]	cycle started by
[CAB]	cycle aborted by
[LRB]	load retrieved by
[LRT]	time user retrieved load ("all data" viewing only) (hh:mm:ss, dd\mm\yy)
(pin not displayed)	

Interpreting the output

To indicate a successful cycle, [STS] displays "Cycle Complete" for 134°C N cycles, or "Calibration Complete" for calibration cycles.

If any cycle parameters are incorrect, [STS] displays ***FAIL***, followed by an Advisory Code. Refer to Appendix 4.

For 134°C N cycles, you can check that:

Temperature [TMX] and [TMN] is between 134-137°C

Time between [TSS] and [TSE] is >3min.20secs

Pressure @134°C = 303.5kPa, @137°C = 331.1kPa

Eschmann SES 2010			
[SRN]	LENA9I1001		
[CYC]	134°N	00069	
[STD]	11:34:17.9	18\10\19	
	Time	kPa	°C
[TSS]	11:41:26.2	318	135.0
	11:42:26.2	321	135.4
	11:43:26.2	321	135.4
	11:44:26.2	322	135.6
[TSE]	11:44:46.2	318	135.4
[TMX]			135.7
[TMN]			135.0
[DCB]	00		
[PTD]	11:45:46.5	18\10\19	
[LRB]	01		
[ETD]	11:47:29.1	18\10\19	
[STS]	Cycle Complete		



14.5 Connecting a Printer

Eschmann can supply an optional printer.

The Printer plugs into a standard domestic mains socket. It connects to the autoclave via an RS232 cable connected to the RS232 port (Fig. 1 item 6).

It is permanently attached to the autoclave, and will print a cycle record as each cycle progresses.

Note: *In order for the printer to operate it must first be activated by an Eschmann or Eschmann trained engineer.*

With careful storage and avoiding sunlight, printed records can be stored for many years.

To order:

REF 8727390 Thermal Printer kit

REF 8727078 Printer Rolls (5)

Appendix 1 : Daily/weekly testing

Manufacturer's testing must be carried out as outlined below.

We recommend a daily Automatic Control Test, and the completion of a Daily/Weekly Log Book (available from Eschmann - 119857 Log Book (N Type)).

All tests must be completed successfully.

Daily Testing

Automatic Control Test (with TST Indicator)

1. Load the autoclave and the TST Indicator as per the manufacturer's instructions.
2. Run the 134°C N type cycle.
3. Remove the TST Indicator at the end of the cycle.
4. Compare the TST Indicator with the Device's instructions.



WARNING!

If errors are detected when running test cycles, or if the TST Indicator indicates a failure, the cause of the failure MUST be corrected and the test must be repeated successfully before the autoclave is used to process loads.

CAUTION

The TST Indicator must only be used and stored as detailed in the Instructions For Use supplied with the Device. Failure to do so could lead to failure of the Device and dangerously misleading results.

This device satisfies the requirement to perform a daily steam penetration test as given in:

HTM01-05 Decontamination in primary care dental practices.

EN 17665-1 Sterilisation of health care products. Moist heat. Requirements for the development, validation and routine control of a sterilisation process for medical devices.

DB 2002(06) Device Bulletin, Benchtop Steam Sterilisers - Guidance on purchase, Operation and Maintenance, Medical Devices Agency, clause 5.3.1.

HTM01-01 (HTM2010/DB2002(06)) testing

If required, further Daily, Weekly, Quarterly and Annual Testing may be necessary to satisfy this guidance. Contact Eschmann for details on the automatic control test.



Fig. 13. TST Indicators (Class 6)

To order:

REF 8793547 TST Indicators (Class 6) - 250 uses

Routine processing

A TST Indicator is to be placed into the centre of the load for every cycle processing instruments.

Appendix 2 : Loading the autoclave

⚠ WARNING!

When loading the autoclave take care not to damage the door or front face of the chamber, especially the door seal and mating face. Damage to these parts can adversely affect performance.

Do not reprocess single-use items.

- Thorough cleaning of all items prior to sterilisation is imperative to ensure instruments can be effectively sterilised. If visible debris or bioburden is not removed prior to sterilisation it will interfere with microbial deactivation and will compromise the sterilisation process.
- Eschmann always recommend the use of a validated, automated cleaning process, i.e. thermal washer disinfectant. Always follow the instrument manufacturer's instructions.
- Ensure all instruments are rinsed thoroughly prior to sterilisation to remove any traces of chemical detergents/disinfectants as these may cause corrosion during the sterilisation process. Ensure instruments are dried prior to being placed in the autoclave.
- Open and unlock all items fitted with hinges or ratchets.
- Dismantle, or only loosely assemble, any multiple-part items to allow steam penetration.
- Once loaded, start the autoclave cycle. Do not leave items awaiting processing in the autoclave, unless set for delayed start (Refer to Section 12.1).
- Check the suitability of non-metallic items for processing in an autoclave. Check that the temperatures are acceptable.
- Load types and load arrangements should conform to that specified below and detailed in Appendix 5 : "Maximum Loads". If in doubt about loads and processing details, contact Eschmann.

Unwrapped solid items only

- The autoclave is ONLY suitable for unwrapped solid items. Refer to Section 1.3 Limitations of Use.
- Place unwrapped solid items loose on trays or into hands-free instrument transfer (HFIT) cassettes.
- DO NOT overload trays or HFIT cassettes (refer to Appendix 5: "Technical data" for details of maximum tray loads and maximum load).
- Avoid bunching items together and ensure all items are positioned so that they do not touch and can drain freely.
- Only use trays or HFIT cassettes supplied by Eschmann.
- Place loaded trays and HFIT cassettes in the sterilising chamber.
- Position all items so that they drain freely and do not trap rising air bubbles.
- Sterilised items are intended for immediate use.

Using the tray lifter

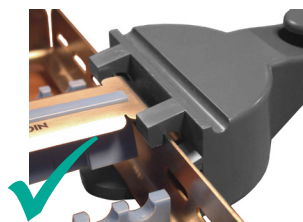
CAUTION

Take care when unloading.
The load and the chamber may be very hot.


Remove the load safely using the tray lifter. Position the half-round base under the item to be lifted and hook on to the item using one of the pairs of top clips:

- Use the upper pair for HFIT (Hand Free instrument Transfer) cassettes - hold at the side or end of the cassette.
- Use the lower pair for trays - not the upper pair.
- Pull the tray carrier to unclip it from the chamber. The tray carrier requires a gentle push to clip it back into the autoclave chamber.

HFIT cassettes




HFIT cassette gripped on its side.

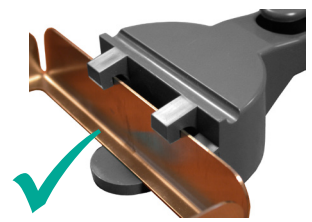


HFIT cassette gripped on its end.

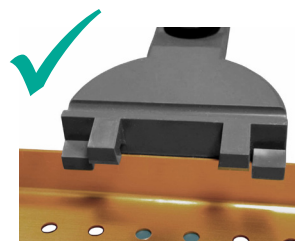Trays



Tray in bottom clips
CORRECT

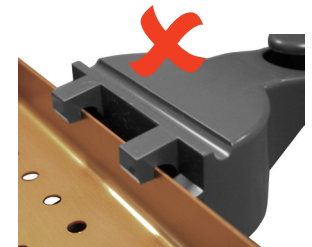

Tray in top clips
INCORRECT

To order:

REF 87-271-21 Tray Lifter

Load capacity


The load should be confined to the usable space.

Maximum load: 4.5kg

*Including the weight of the trays

Typical loads

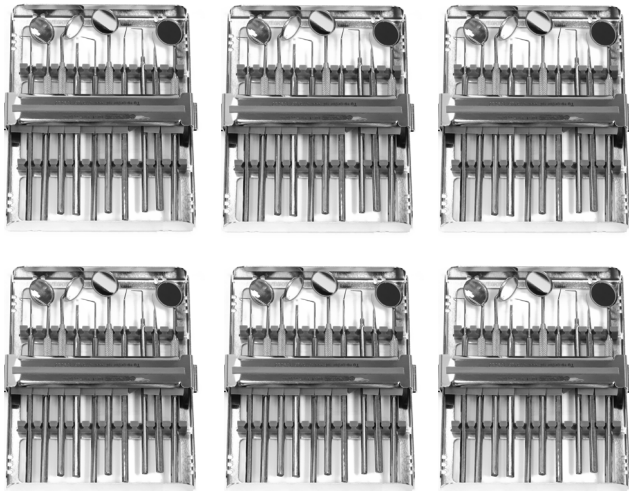
Example 1: 6 x HFiT10 Cassettes (60 Instruments)



Loading configuration:

Top row:	1 x HFiT10 Cassette
Upper middle row:	2 x HFiT10 Cassettes
Lower middle row:	2 x HFiT10 Cassettes
Bottom row:	1 x HFiT10 Cassette

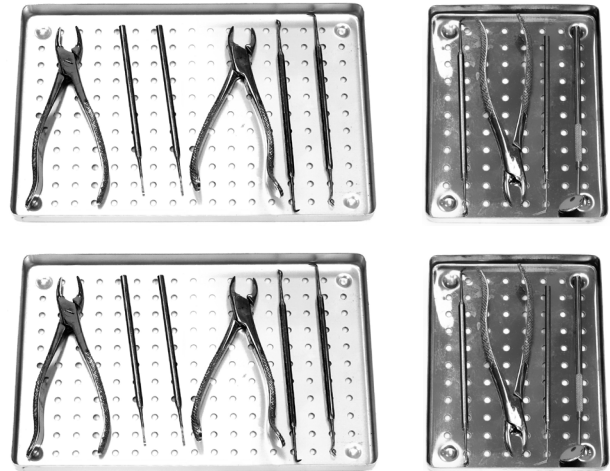Example 2: 2 x Standard Trays and 2 x Half Trays



Loading configuration:

Top row	1 x Half Tray
Upper middle row:	1 x Standard Tray
Lower middle row:	1 x Standard Tray
Bottom row:	1 x Half Tray

Maximum load per Standard Tray: 1.5kg*

Maximum load per Half Tray: 0.75kg*

*Including the weight of the tray

The SES 2010 Autoclave is used to sterilise medical instruments and consideration should be given to hazards inherent in the nature of the materials to be autoclaved. Any consumable (e.g. door seals, trays and other chamber furniture) should follow local biohazard procedures for cleaning and decontamination prior to disposal.

This autoclave is used to decontaminate potentially pathogenic material from instruments used in the Health sector. Eschmann are able to advise on how to recycle and dispose of the product correctly.

Appendix 3 : Maintenance

Cleaning and care

CAUTION

Even in the concentrations found in tap water, chlorine can cause damage to the autoclave.

Disconnect from the mains electrical supply before cleaning the Autoclave.

Do not use flammable liquids, abrasive powders, chemicals, or solutions containing chlorine to clean the autoclave

Note: For guidance on clean steam management contact Eschmann.

- Check the door seal and chamber face for damage every day.
- Clean the door seal and chamber face every day using Eschmann recommended wipes. DO NOT use tap water.
- Keep the chamber and chamber trays clean.
- Clean and sterilise the tray lifter regularly (e.g. once a week).
- Clean the outside of the autoclave using Eschmann's range of recommended cleaning products.

CAUTION

In common with other systems containing static water reservoirs, the water in this autoclave can become contaminated over a period of time and should be treated as a potential risk of infection.

- As a minimum, drain the water reservoir at the end of each day.
- Check the drain tube (Part No. 115036) and the connector regularly. Replace if damaged.
- Always leave the chamber empty and the door ajar when not in use.
- Fully drain the reservoir if preparing the autoclave for transit or storage.
- In-service inspection and testing of electrical equipment must be carried out routinely in accordance with local legislation.

WARNING!

Should the autoclave ever leak, disconnect from the mains electrical supply, drain the reservoir and call Eschmann.

Autoclave fuses

Fuse x 2 (Fig. 1, pg 4).

Refer to Appendix 5: "Technical data".

1. Switch off and disconnect the autoclave from the electrical mains supply.
2. Remove the fuse cover using a screwdriver or small coin (turn anti-clockwise).
3. Inspect and renew the fuse if necessary.
4. Refit the fuse cover (turn clockwise).

Mains cable fuse

Fuse x 1

Refer to Appendix 5: "Technical data".

1. Switch off and disconnect the autoclave from the electrical mains supply.
2. Remove the fuse from the mains plug.
3. Inspect and renew the fuse if necessary.
4. Reconnect the plug.

Waste electrical & electronic equipment (WEEE) regulations

The aim of the WEEE Regulations is to reduce the amount of waste going to landfill.

 All Eschmann products that must be recycled in accordance with the WEEE Regulations are marked with the "wheelie bin" symbol opposite.

What to do

Please contact us when one of our products, marked with the symbol, reaches the end of its working life. We will be able to advise on how to recycle and dispose of the product correctly.

Note: If we request that the product is returned to Eschmann, it must be decontaminated first. We will request a certificate.

Under the WEEE Regulations, manufacturers are held responsible for recycling waste electrical and electronic equipment (WEEE) placed on the market after 13 August 2005 that has reached the end of its working life. The regulations also place obligations to comply on distributors, retailers and end users of the equipment.

Service, Calibration & Certification

Please contact Eschmann to discuss servicing and certification requirements. The autoclave must be serviced and calibrated by Eschmann or Eschmann trained engineers at 12 month intervals, or when the Service Required screen is displayed, whichever occurs sooner. Refer to Section 10.2. The autoclave must be periodically examined and certified against the Pressure Systems Safety Regulations 2000 (PSSR). Eschmann recommends this is performed 18 months from manufacture and every 14 months thereafter.

Appendix 4: Trouble shooting

Trouble shooting table

Fault	Action
Autoclave will not power-up	Check the mains supply, mains connections, autoclave fuses and mains cable fuse. If fault cannot be found, contact Eschmann.
Autoclave powers-up but the screen is blank	Power autoclave off. Remove any item plugged into the USB port (Fig. 1, pg 4). Power autoclave on. If unsuccessful, contact Eschmann.

Advisory code table

The autoclave continually checks that it is able to run a cycle and sterilise the load to the required standard.

If it encounters a problem it will inform you via the Advisory Screens:

- Water Management Screens.
- Service Screens.
- Advisory Screens:
 - POWER FAILURE
 - USER STOPPED CYCLE
 - SYSTEM STOPPED
 - ERROR

Refer to Section 11 Advisory Screens

All Advisory Screens are identified internally using an Advisory Code. These are shown in the Advisory Screens Table below.

Most of the problems referred to by the Advisory Screens can be rectified, allowing the cycle to continue.

In most cases, **press OK to clear the Advisory Screen** and follow the advice in the table below.

Note: *The advisory screens that cannot be cleared by the operator are the ERROR screens. Contact Eschmann. Eschmann require the Advisory Code from these screens to determine the problem. See end of table 1.*

Table 1: Advisory Codes

Advisory screens	Advisory code	Description	Action
POWER FAILURE	1	Power Failure	Check for intermittent mains supply and mains connections. Run cycle again.
USER STOPPED CYCLE	2	User Abort	Run the required cycle.
SYSTEM STOPPED	3	Door Lock Fail	Check for obstructions. Run cycle again.
	11	Temperature Synchronisation Fail	Power autoclave off / on. Run cycle again.
	15	Arbiters Disagree	Power autoclave off / on. Run cycle again.
	16	Invalid Cycle Type	Power autoclave off / on. Run cycle again.
	19	Temperature and Pressure Disagree	Visually check door seal. Clean mating face. Run cycle again.
	20	Temperature Sensors Disagree	Power autoclave off / on. Check chamber and load is dry. Run cycle again.
	21	T1 High During Sterilisation	Power autoclave off / on. Run cycle again.
	23	T1 Low During Sterilisation	Power autoclave off / on. Run cycle again.
	25	Not Enough Water in Reservoir	Fill with fresh water. Run cycle again.
	26	Too Much Water in Reservoir	Empty waste water. Run cycle again.
	33	Door Not Closed	Power autoclave off / on. Run cycle again.
	34	Door Lock Fail	Power autoclave off / on. Run cycle again.
	36	PR and/or GI Not In Running	Power autoclave off / on. Run cycle again.
	41	Control System Error Occurred	Power autoclave off / on. Run cycle again.
42	Error Wrong Hold Time	Power autoclave off / on. Run cycle again.	

Advisory screens	Advisory code	Description	Action
	43	Missed First Process Assessment Number	Power autoclave off / on. Run cycle again.
	44	Process Assessment Sequence Error	Power autoclave off / on. Run cycle again.
	45	Error Unexpected End of Process evaluation List	Power autoclave off / on. Run cycle again.
	48	Error Process Assessment Failed	Power autoclave off / on. Run cycle again.
	49	Drying Failed	Power autoclave off / on, run cycle again. Run without drying.
	155	Communications Failure	Power autoclave off / on. Run cycle again.
	156	Communications Request Failure	Power autoclave off / on first, then press OK to clear the SYSTEM STOPPED screen. Run cycle again.
	157	Communications Remove Failure	Power autoclave off / on first, then press OK to clear the SYSTEM STOPPED screen. Run cycle again.
	160	System Synchronisation Failure	Power autoclave off / on. Run cycle again.
	175	Pressure Not Equalised	Power autoclave off / on. Run cycle again.
	176	Safety Relay Fail	Visually check door seal. Clean mating face. Run cycle again or Leave autoclave to cool. Safety relay resets automatically. Run cycle again.
	201	Enclosure Temperature Too High	Leave autoclave to cool. Run cycle again.
	202	Process Assessment Temperature Too High	Power autoclave off / on. Run cycle again.
	203	Process Assessment Temperature Too Low	Power autoclave off / on. Run cycle again.
	204	Process Assessment Pressure Too High	Power autoclave off / on. Run cycle again.
	205	Process Assessment Pressure Too Low	Power autoclave off / on. Run cycle again.
	209	Process Assessment Exceeded Time Limit	Power autoclave off / on. Run cycle again.
	210	Process Assessment Sterilisation Temperature Too Low	Power autoclave off / on. Run cycle again.
	211	Process Assessment Sterilisation Temperature Too High	Power autoclave off / on. Run cycle again.
	217	Data Corrupt	Power autoclave off / on first, then press OK to clear the SYSTEM STOPPED screen. Run cycle again.
	223	Touch Screen Continuously Pressed	Power autoclave off / on. If persistent contact Eschmann.
	226	Cycle Check Failed	Power autoclave off / on. Run cycle again.

Advisory screens	Advisory code	Description	Action
	228	Start Up Test Failed	Power autoclave off / on first, then press OK to clear the SYSTEM STOPPED screen. Run cycle again.
	231	Calculated Temperature Too High	Power autoclave off / on. Run cycle again.
	232	Calculated Temperature Too Low	Power autoclave off / on. Run cycle again.
	233	Difference between Calculated Temperature and Chamber Temperature too large	Power autoclave off / on. Run cycle again.
	909	Chamber Temperature Sensor	Power autoclave off / on. Run cycle again.
	910	Chamber Pressure Sensor Fault	Power autoclave off / on. Run cycle again.
	913	Enclosure Temperature Sensor Fault	Power autoclave off / on. Run cycle again.
	916	Unknown Sensor Fault	Power autoclave off / on. Run cycle again.
	919	Wrong unit type	Press OK to clear the system stop screen. If persistent contact Eschmann.
	920	Unknown unit type	Press OK to clear the system stop screen. If persistent contact Eschmann.
	921	CO not fully configured	Press OK to clear the system stop screen. If persistent contact Eschmann.
	996	Water quality system fault	Press OK to clear the system stop screen. If persistent contact Eschmann.
	997	Protective system temperature fault	Power autoclave off / on. Run cycle again. If persistent contact Eschmann.
	998	Protective system temperature fault	Power autoclave off / on. Run cycle again. If persistent contact Eschmann.
ERROR	18	Protective system pressure fault	Power autoclave off / on. Run cycle again. If persistent contact Eschmann.
	177	Failed To Update Calibration Record	Contact Eschmann.
	178	Calibration Invalid	Contact Eschmann.
	206	Safety Valve Failure	Contact Eschmann.
	208	Safety Valve Failure	Contact Eschmann.

Appendix 5 : Technical data

POWER SUPPLY			
Supply Voltage	Supply 230V (±10%) or 220V (-6%, +15%) See rating plate. 50/60 Hz a.c. ~ For use with alternating current		
Fuse Rating (panel)	2 x F10A, 250V (Part No. MUW218616)		
Fuse Rating (plug)	13A		
SAFETY CONDITIONS			
Electric Shock Protection	Class 1		
Pressure Relief Valve	Set pressure: 2.85 barg		
STERILISING AND TEST DATA			
Sterilising Time	3 minutes 20 seconds at 134/137°C The sterilisation time can be increased to a maximum of 30 minutes by an Eschmann trained engineer.		
Drying Time (approximate)	17 minutes		
Total cycle time	11 mins (Empty, no drying)		
Operating Pressure	303.5 kPa (abs) at 134°C, 331.1 kPa (abs) at 137°C		
MAXIMUM LOADS			
Maximum load per standard tray	1.5 kg		
Maximum weight per half tray	0.75 kg		
Maximum total load	4.5 kg		
WATER AND RESERVOIR			
Water Reservoir capacity/ Use per cycle	2500 ml/ 600ml (Proportion re-cycled)		
Temperature of drained water/ flow rate	81°C/ 1.2 l/min		
Heat transmitted by autoclave in one hour	720 kj		
DIMENSIONS AND WEIGHT			
Depth	468 mm		
Width	460 mm		
Height	360 mm		
Usable chamber space	190 x 110 x 288mm		
Weight (empty)	26.5kg		
Weight (full reservoirs and full load)	33.5kg		
SHIPPING DIMENSIONS AND WEIGHT (includes unit, accessories and packing)			
Depth	605mm		
Width	655mm		
Height	468mm		
Weight	31.5kg		
STANDARDS COMPLIANCE			
BS EN ISO 15223-1	BS EN 13060	BS EN 61326	BS EN 13445
BS EN ISO 20417	BS EN ISO 14971	BS EN 62366	BS EN 61010
ENVIRONMENTAL REQUIREMENTS (contact Eschmann if required to operate outside these conditions)			
Transport and Storage: (Note: Water must be drained)			
Ambient temperature range	-10°C to +40°C		
Relative humidity range	30-70%RH (Non-condensing)		
Atmospheric pressure range	700 millibars to 1060 millibars		
Operation:			
Location	Indoor use only		
Ambient temperature range	+5°C to +40°C		

Relative humidity	Maximum relative humidity 80% up to 31°C decreasing linearly to 50% relative humidity at 40°C
Altitude	Maximum 2000m (800 millibars)
MISCELLANEOUS	
IFU applicable to serial numbers from:	LENA9J0202
Basic UDI-DI (BUDI-DI) Identifier	506069061AUTOCLAVEF9
Technical Lifetime	7 years when used in accordance with these instructions, properly maintained and repaired.
Servicing, calibration and testing interval	Annual Eschmann products must be serviced by Eschmann, or Eschmann trained engineers only . Failure to do so may invalidate the warranty.
Door Seal Replacement Frequency	500 cycles
Autoclave Sound Level	58 dBA
Sounds	Autoclave Ready: 1 beep, 300 ms Keypress: 1 beep, 100 ms Cycle Start: 1 beep, 100 ms Cycle End: 1 beep, 500 ms Cycle Failed: 5 beeps, 100 ms Incorrect PIN: 2 beeps 100 ms
Chamber Label	Chamber volume: 10.6 litres Max design pressure: Pd Max 2.85 bar Min design pressure: Pd Min -1.00 bar Design Temperature: Td 140°C Max Allowable Pressure: Ps 2.85 bar Safety Valve Set Pressure: 2.85 bar Hydrostatic Test Pressure; 4.15 bar
Pollution degree	2

SAFETY FEATURES

Eschmann autoclaves are designed to be safe and effective. No changes should be made, except by an Eschmann Engineer. In particular, the following safety features must not be interfered with, circumvented or overridden:-

- Door closed (position) sensor and door lock (position) sensor. These prevent a cycle starting if the door is not closed and locked.
- Pressure Relief Valve (safety valve) set at 2.85+10% barG.
- Independent cycle control and monitoring. There are two microcomputers independently monitoring each other and other autoclave components. Either is able to stop a cycle and place the autoclave into a safe condition in the event of a failure. In the case of a venting valve failure the autoclave may not be able to eject water and steam (i.e. pressure will remain until cool) and the door will remain locked until the pressure / temperature has dropped to a safe level.
- Overheat Control. The microcomputer operates in conjunction with an independent thermostat to protect the heating element from overheating.
- Pressure door switch: Prevents the door being opened when the residual chamber pressure is too high by removing power to the door interlock so that it locks the door closed.

BS EN 13060 Type tests


In accordance with BS EN 13060 the recommended type tests for N type cycles are:

Type Tests	
Dynamic steriliser chamber pressure	
Air Leakage	
Empty chamber	✓
Solid load	✓
Small porous loads	
Full porous load	
Narrow Lumen	
Dryness, solid load	
Dryness, porous load	
Non-condensable gases	✓

Load chart

	Cycle Data
Sterilisation Temperature	134-137°C
Sterilisation Pressure	303.4-331.1kPa (abs)
Duration of plateau phase	3min.20secs
Duration of drying phase	17 minutes
Total Cycle Duration (excludes drying time)	Empty: 11 min Full load: 14 min (Full load, cold machine: 21min)
LOAD TYPE :	Solid (Unwrapped) loads only
MAXIMUM LOAD : Weight Limits	4.5kg

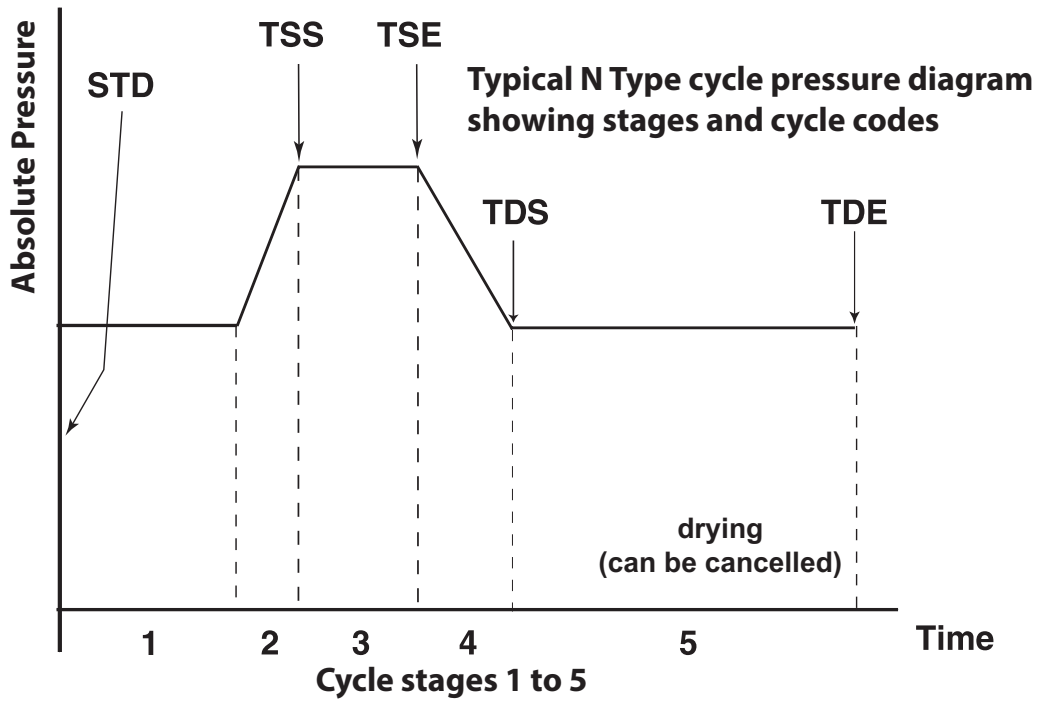Detailed N type cycle information



Progress Indicator	Cycle stage	Stage Detail
	1 & 2	(1) Filling and (2) heating to sterilisation temperature.
	3	Sterilising, temperature and pressure held: 134-137°C, 303.5 - 331.1kPa
	4	Discharge phase.
	5	Drying phase (Selectable, default is no drying)

Appendix 6 : Accessories

For prices and ordering call Eschmann, telephone 01903 753322

Visit Eschmann.co.uk for details of the full range of decontamination equipment, accessories and consumables.

Autolog Wireless Cycle Data Manager - for more advanced automated cycle logging and digital generation of daily/weekly testing, in real time, upgrade to the Autolog

8701028



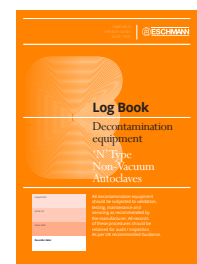
Essential Test Devices

8793547 TST Indicators - Class 6 (box of 250)



Logbook

119857



RO Water System

8785154UVCSHW



Standard Trays (28cm x 18cm)

8724451	Aluminium Silver	8724454	Aluminium Green
8724452	Aluminium Blue	8724455	Aluminium Yellow
8724453	Aluminium Red	8724456	Aluminium Turquoise
8725450	Stainless Steel		

Examination/ Half Trays (18cm x 14cm)

8724461	Aluminium Silver	8724464	Aluminium Green
8724462	Aluminium Blue	8724465	Aluminium Yellow
8724463	Aluminium Red	8724466	Aluminium Turquoise
8725460	Stainless Steel		




HFIT Cassettes

HFIT 10 'Examination' Stainless Steel Cassette (13.5cm x 18cm)

8722812	Stainless Steel Blue	8722813	Stainless Steel Red
8722814	Stainless Steel Green	8722815	Stainless Steel Yellow
8722816	Stainless Steel Silver	8722818	Stainless Steel White

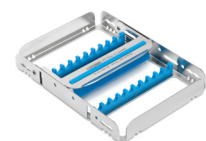HFIT 5 & 20 sizes also available



Spares and Accessories

8727121	Tray Lifter	MUW216747	Door Seal Kit
115036	Spare Drain Tube	115664	USB Memory Stick
MUW218616	Mains Fuse		

Service

The Autoclave must only be serviced by an Eschmann or Eschmann trained engineer. We also recommend a program of ongoing service and maintenance using only Eschmann spare parts.

Eschmann products are supported by a worldwide network of fully trained engineers, offering high quality Eschmann spare parts. For further information on the range of Service Contracts available, please contact your local Eschmann representative.

Eschmann can be contacted during normal office hours. Please quote the model and serial number exactly as printed on the product label. Please use the contact details below.

IMPORTANT: The design of the autoclave pressure vessel is certified by a third party accredited test house to International Standards. In order to ensure safety and to comply with UK and/or International regulations, the vessel and fittings should be inspected by a competent person at regular intervals. We recommend that this is carried out at least once every 14 months by an Eschmann or Eschmann trained engineer.

Safety Reporting

Any serious incident that occurs with this medical device should be reported to Eschmann Technologies Limited and your local Competent Authority

Warranty

This product and purchased accessories are warranted for a minimum period of 12 months to be free from defects in materials and workmanship at the time of delivery.

Eschmann will be under no liability for any defect arising from fair wear and tear, negligence, wilful damage, misuse, abnormal working conditions, failure to follow the manufacturer's instructions, unauthorised alteration or repair of hardware, unauthorised or accidental alteration of software or configuration, lost profits, commercial loss, economic loss, or loss arising from personal injury. We may, at our discretion, raise a charge for any faults repaired that fall outside the warranty cover. Where charges are necessary, replacement parts will be charged at manufacturers' list prices and labour will be charged at the prevailing hourly rate. Repairs performed by Eschmann carry a 3-month parts and labour warranty.

The details in this manual are correct at time of going to press and supersede details previously published elsewhere. Information, technical specifications and statements made in this publication may be subject to change without prior notice

EC Authorised Representative

Casus Europe B.V. Lange Viestraat 2b, 3511 BK Utrecht,
The Netherlands

Eschmann Technologies Limited

Eschmann House, 15 Peter Road, Lancing,
West Sussex, BN15 8TJ United Kingdom

t: +44 (0)1903 753322

e: info@eschmann.co.uk

w: eschmann.co.uk

Ingenuity in
infection control
since 1830

